
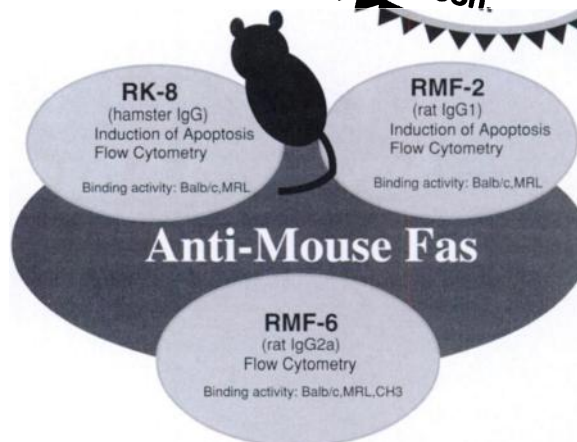


Anti-Mouse Fas & Anti-Human Fas

Monoclonal Antibodies

MBL

ELISA Kit
for measuring
soluble Fas
will be available
soon.

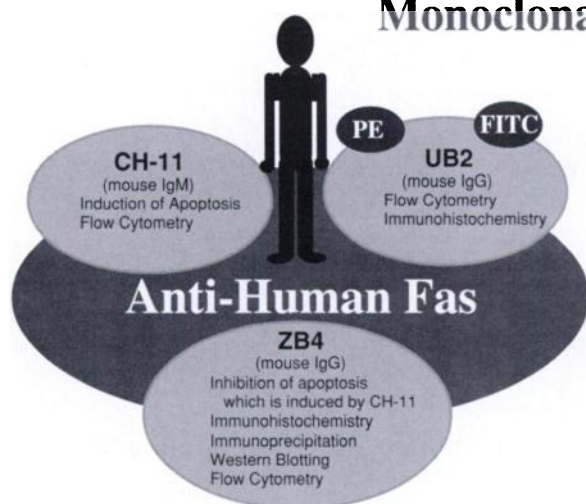



	clone		Induction of apoptosis	Inhibition of apoptosis	Flow Cytometry	Immunohistochemistry	IP	WB
	CH-11	purified	●		●			
	UB2	purified			●	●		
		FITC conj.			●			
		PE conj.			●			
	ZB4	purified		●	●	●	●	●
	RK-8	purified	●		●			
	RMF-2	purified	●		●			
	RMF-6	purified			●			

MEBSTAIN Apoptosis kit


Principle

In cells in which apoptosis occurs, the chromatin DNA is cut by endonuclease at linker DNA site between nucleosomes. Then, DNA fragments, which are a number of multimers of nucleosomal units exist in nuclei.

In TUNEL method, 3'-OH DNA ends generated by DNA fragmentation is nick end labelled with biotinylated dUTP, introduced by terminal deoxytransferase, and then stained using avidin-conjugated FITC to allow specific staining.

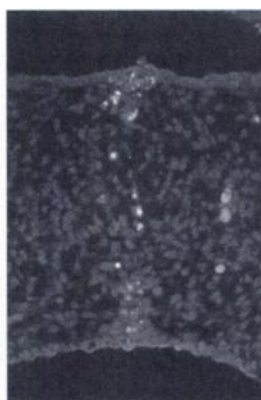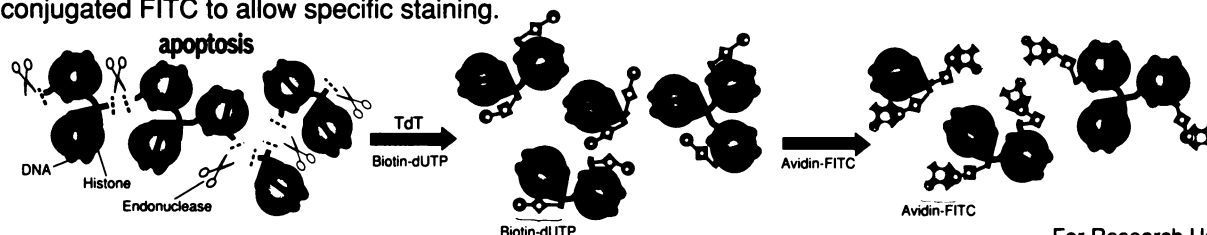

Fetal mouse midpalate fusion

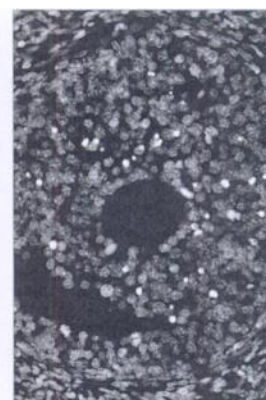

Normal ovary (mouse)

Photos was offered by
Dr. Chisato Mori, Dept. of Anatomy,
Faculty of Medicine, Kyoto University.

For Research Use Only



MBL International Corporation
440 Arsenal Street, Watertown, MA 02172, USA
TEL: (617) 926-6964 FAX: (617) 926-6954



MEDICAL & BIOLOGICAL LABORATORIES CO., LTD.
5-10 Marunouchi 3 chome, Naka-ku, Nagoya 460, JAPAN
TEL: (052) 971-2081 FAX: (052) 971-2337



Essentials for Immunohistochemistry

Counterstains

Hematoxylin — Based on Gill's formula, our hematoxylin is a progressive stain, free of mercury or alcohol and especially designed for immunohistochemical applications.

Cat. No. H-3401 • 500 ml

Methyl Green — Our improved formula provides nuclear staining with a simple protocol. This light green counterstain is ideal for black and white photography and for single or multiple labeling procedures.

Cat. No. H-3402 • 500 ml

Nuclear Fast Red — An easy to use counterstain for blue, brown or black enzyme reaction products.

Cat. No. H-3403 • 500 ml

Antigen Unmasking Solution

This citric acid-based formula was developed to reveal antigens in formalin-fixed, paraffin-embedded tissue sections using high temperature treatment procedures. Provided as 100x stock.

Cat. No. H-3300 • 250 ml

VECTABOND™ Reagent

This tissue section adhesive has become the standard for use in bonding tissue sections to glass slides. Employed for routine paraffin or frozen sections or for antigen unmasking or *in situ* hybridization procedures using high temperature, this reagent avoids problems associated with glue or denatured protein adhesives. Provided as a 50x stock sufficient for treating about 500 slides.

Cat. No. SP-1800 • 7 ml

VECTOR LABORATORIES, INC.

30 Ingold Road
Burlingame CA 94010
(415) 697-3600 Fax: (415) 697-0339

VECTOR LABORATORIES, LTD.

16 Wulfric Square, Bretton,
Peterborough PE3 8RF U.K.
(01733) 265530 Fax: (01733) 263048

plus over 25 distributors worldwide